
Dranspose

A constrained map-reduce live analysis pipeline for fast experimental feedback


Felix Engelmann, Scientific Data fe-research@nlogn.org

25th September 2024

Experimental Feedback Loop

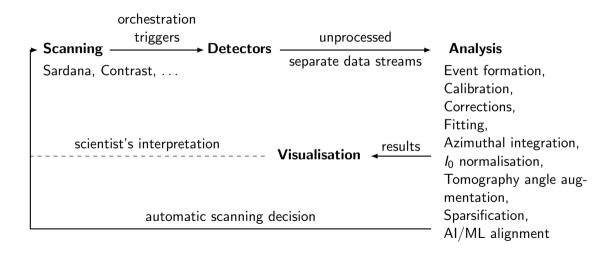

Most Existing Infrastructure

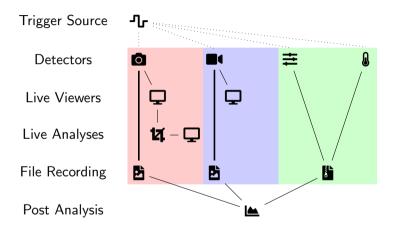

Limitations of Live Analysis

Single Data Stream

- only one detector data available
- ▶ simple tools, e.g. azint, crop, time integration



Beyond Limits

normalisation to I_0 , sorting by motor position

Additional Data Sources

Devices

PandAbox



XSpress3



Basler




AlbaEM


Software

- Sardana
 - Icepap
 - ► Meta-data
 - Filename
 - Snapshots
- Contrast
- Bluesky


Event Formation

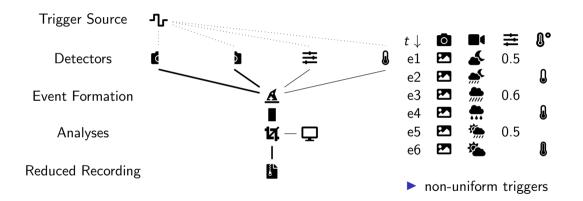

Event Formation Bottleneck

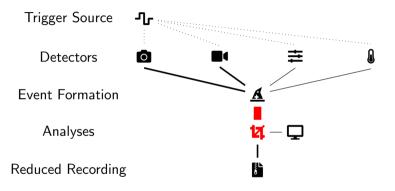

Event Formation Bottleneck

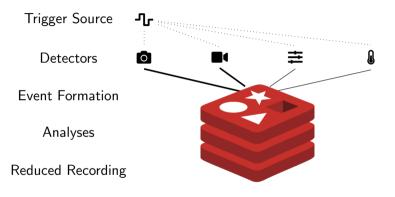

Existing Solutions

Apache Flink

- exactly-once guarantee
- many small events

Hummingbird

- XFEL development
- tailored to pulsed facilities


OnDA / HiDRA

- DESY developments
- maintained?

EWOKS

- large ecosytem
- data ingress from blissdata

Parallelisation by Event, not Stream (map)

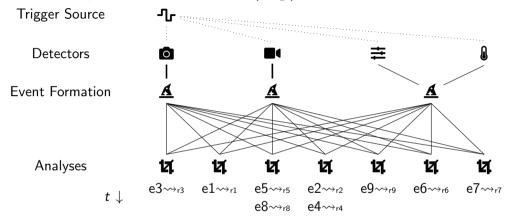

Balancing Constraints

- event 2n and 2n + 1 to same worker
- event m to worker debug


Missing Inter Event Analysis

- time integration
- long temporal correlations


Sequential Reduce

Operations at Acquisition Speed

- append to list
- sum

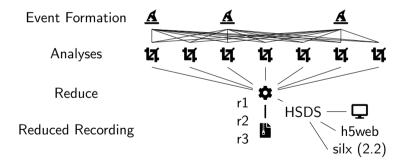


Behind the Scenes

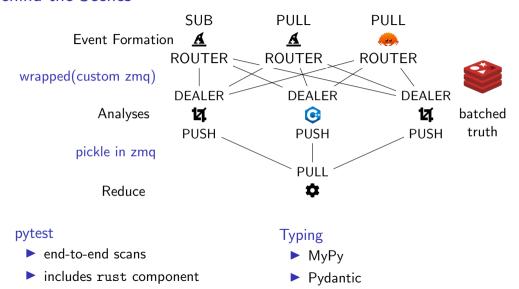

Performance Estimates

K8s cluster of 10: dual Xeon® Gold 6326, 100 Gbit/s Mellanox SR-IOV, 465 GB RAM

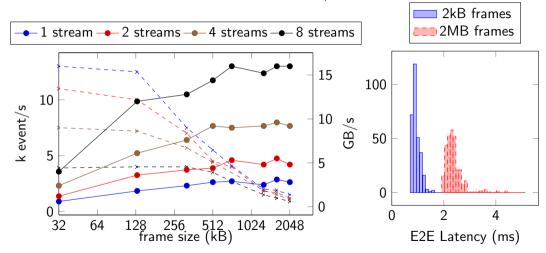

limits on around 2 GB/s per stream or 12 kHz at a few ms

Development: Testing

Recording

- ingesters optionally write all streams to disk (sequential cbor dumps)
- common experimental setup

Replay

- from recorded stream dumps
- from hdf5 files

Local

run file-based ingesters, workers and reducer locally


Case Study: NanoMAX

Sources

- ► XSpress3
- Contrast (Motor Positions)

Analysis / Output

- ► PyMCA fluorescence fitting
- Concentration Map

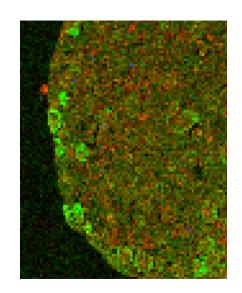

Case Study: DanMAX

Sardana Basler Pilatus Orca PandAbox

Beam CoG, Beam Alignment

Pileup File

Encoder angle



python -m pip install dranspose docker-compose, helm

https://dranspo.se

or come to Max IV, we are hiring: paul.bell@maxiv.lu.se